

# Hotelo Dental X-ray film & "DAYLIGHT" Processing System



*You can diagnose  
by only 30 sec or more.*

## "INSTANT FILM"

(ISO speed D)



DIF size 2 30×40mm  
DIC size 1A 24×30mm  
DIK size 4A 54×70mm

Film packet  
with a slot  
which guides easily to  
insert the nozzle of  
"PUSHER"

## "DQD" solution



Monobath  
solution for rapid  
developing/fixing  
of  
"INSTANT FILM"

## "PUSHER"



Fixed volume injector  
with a nozzle  
attachable  
to the top of  
"DQD" bottle



**HANSHIN TECHNICAL LABORATORY, LTD.**

4-18, Kubocho, Nishinomiya, 662-0927 Japan

PHONE: 81798 336321 FAX: 81798 358113

Quality system of dental X-ray film is certified by TÜV CERT.





# "DAYLIGHT" Processing system for Speed-D INSTANT FILM

## INSTANT FILM

High contrast Speed-D film gives a plenty of information for diagnosis.

Identification film number from 1 to 100 appears on all developed films.

Suitable too for normal developing in automatic processor or in dark-room.


## DQD SOLUTION

Monobath solution for in-packet processing 30 sec at 25°C.

Because of latent high power of the solution, it is able to shorten exposure time to speed-E level in a case of processing at higher temperature (~35°C).


## DEVICES

"PUSHER" and relative devices available in option to support simple process.

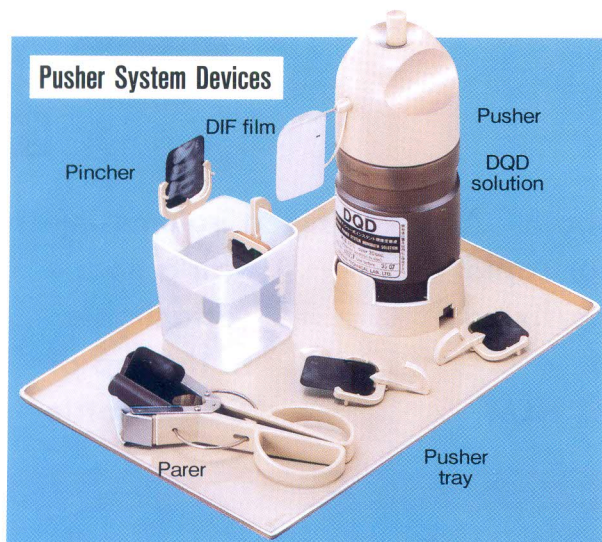

## PROCESSING

"PUSHER" system gives high-quality images through simple operation.

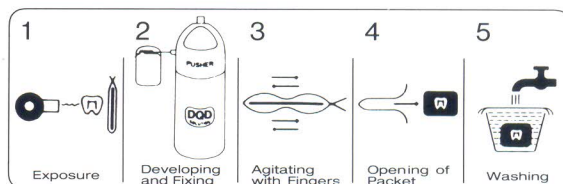

**HANSHIN TECHNICAL LABORATORY, LTD.**

4-18, Kubocho, Nishinomiya 662-0927, Japan.

Distributor:

**N. K. LUCK SDN. BHD.**

Co. No. 75176-V  
91A, JALAN PUDU LAMA  
50200 KUALA LUMPUR  
MALAYSIA

TEL: 03-2070 5622

FAX: 03-2072 3006

E-MAIL: nk Luckkl@tm.net.my